



# Orbital wall prosthesis Codubix® ORBITAL WALL, Codubix® ORBITAL WALL 3D, Codubix® ORBITAL WALL 3D adapted

Prosthesis of orbital wall bone are knitted of polypropylene yarn. High endurance, low specific weight, hydrophobic properties contribute to the lack of liquid and physiological fluid absorption, non-toxicity and chemical inactivity make the prosthesis easy to be identified with the natural bone as for the physical properties. The prostheses can be easily modelled during implantation.

Application:



- for filling orbital wall defects (zygomatic-maxillary-orbital fractures, isolated fractures of orbital wall) resulting from fractures of the central facial part of the skull after injuries and for other reasons

Types:

- sterile (steam sterilised)
- standard version
- 3D version - right and left - anatomical shape
- adapted version

Size and packing:

- standard version: 40 mm x 50 mm
- 3D version: left and right in the sizes S, M, L
- packed in a double gusseted paper-film pouch and then in a box by 1 pc

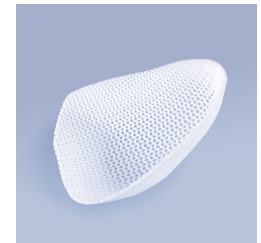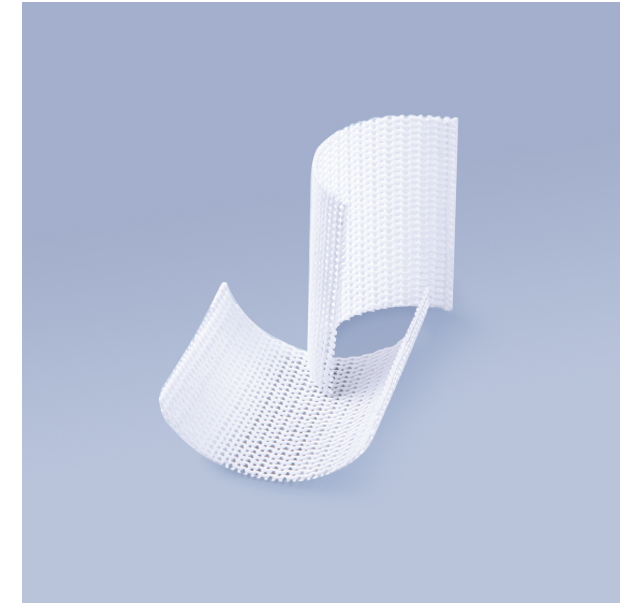

**MD**

Sterile devices are certified by the Notified Body PCBC SA and bear the CE mark with number 1434

>>