Matoset dialysis connection set – application way

Example of a dialysis connection set MA-991-ZEST-487.

1




Before starting the procedure, wash your hands (the specialist and the patient).

2



To perform the procedure, prepare a disinfectant, a tourniquet, suitable cannulas and sets the start and end the procedure.

3

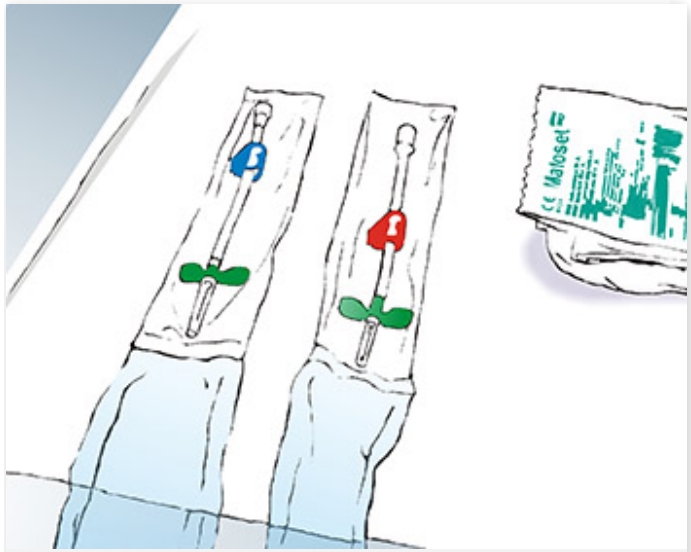

Take out the sterile drape and place it beneath the patient's forearm.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026


Matoset dialysis connection set – application way

4



»

5

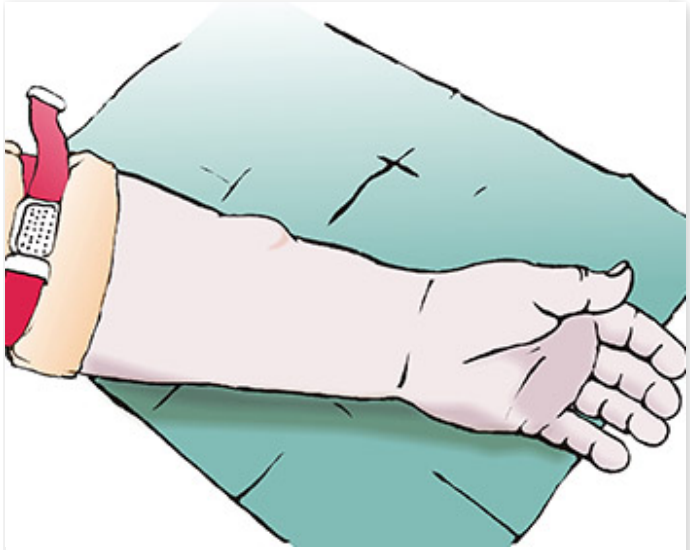

Then put one glove on the left hand, and the other one put on the sterile part of the opened packaging.

6




»

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

Matoset dialysis connection set – application way

7

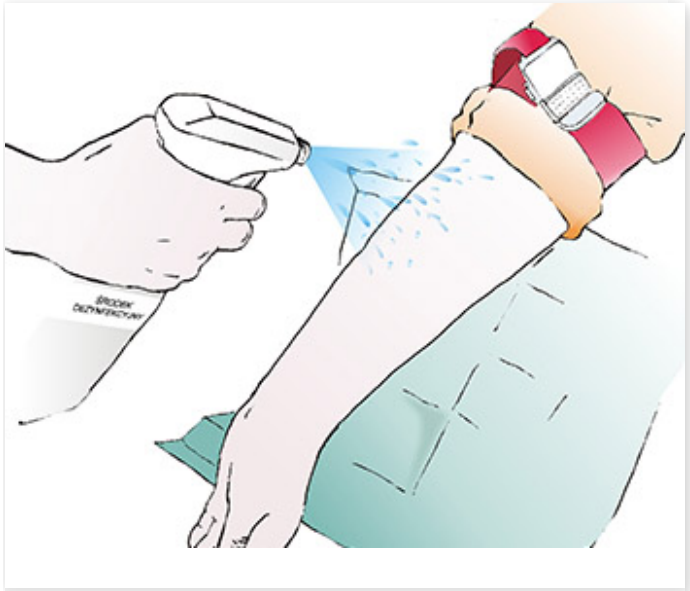

Disinfect the patient's forearm using the prepared disinfection liquid and the sterile swabs from the set.

8



»


9



»

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

10

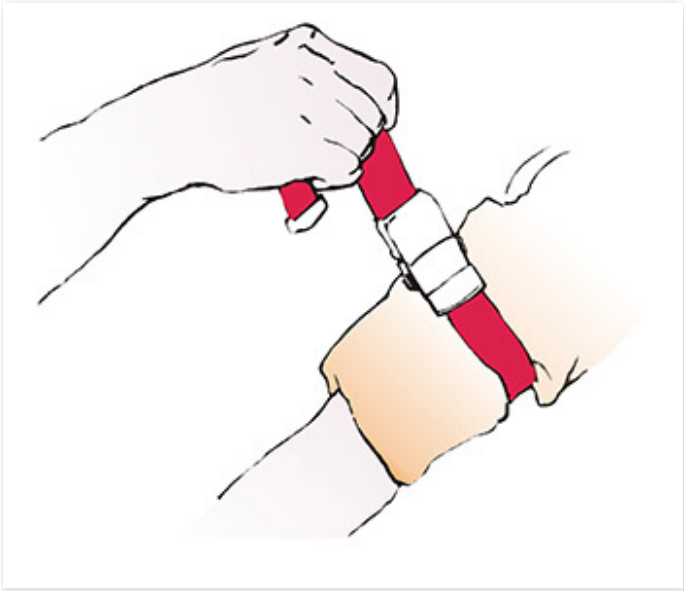


Tighten the tourniquet.

11



Put on the right glove.


12

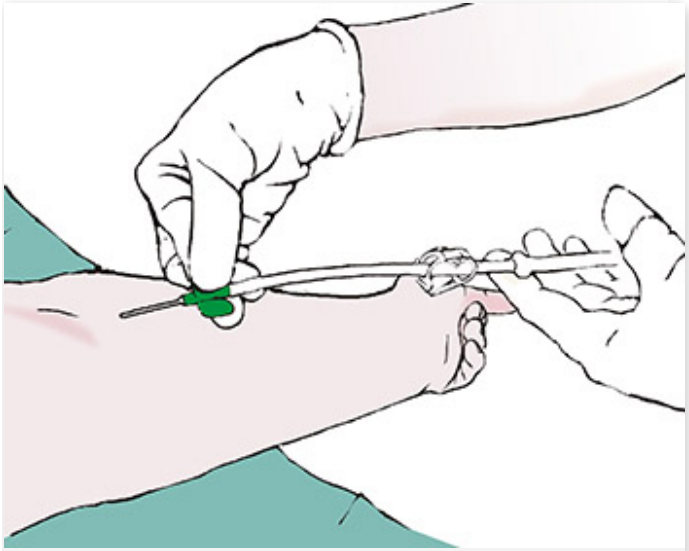

Place the cannulas into the veins.


©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

13

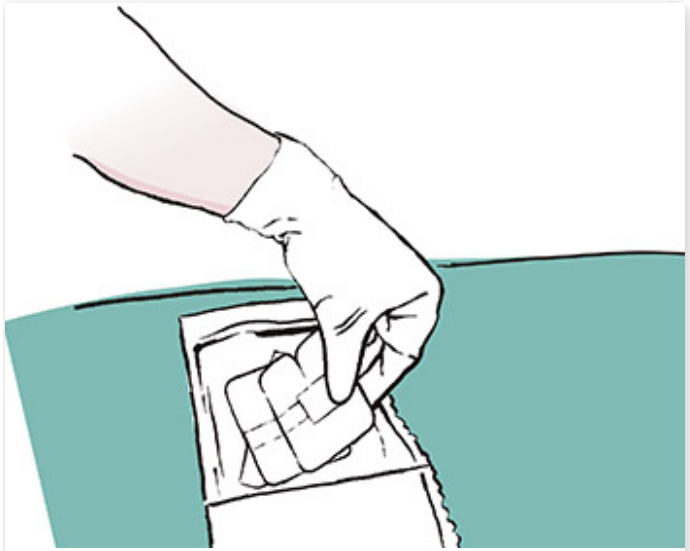


After inserting every cannula fix them with special dressings for fixing cannulas (also in the set).

14

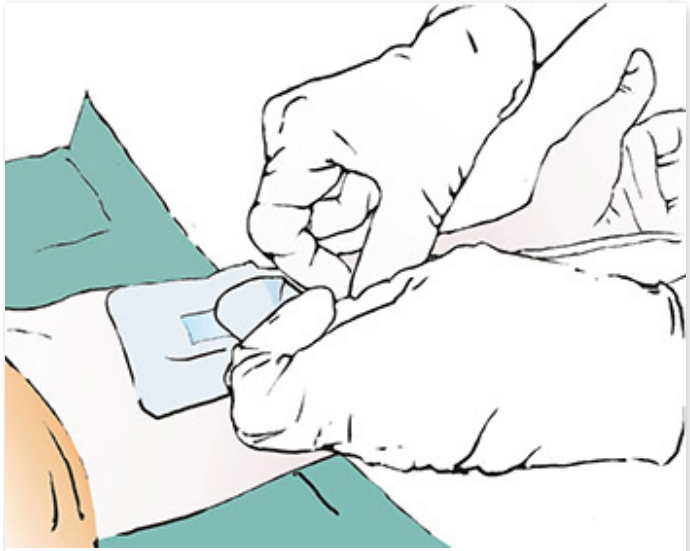

»


15

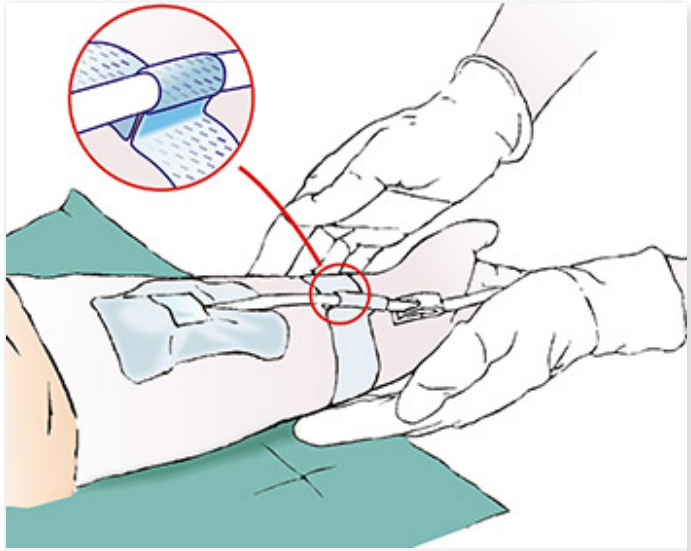

At the end secure cannula drains with the nonwoven adhesive tapes.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

16

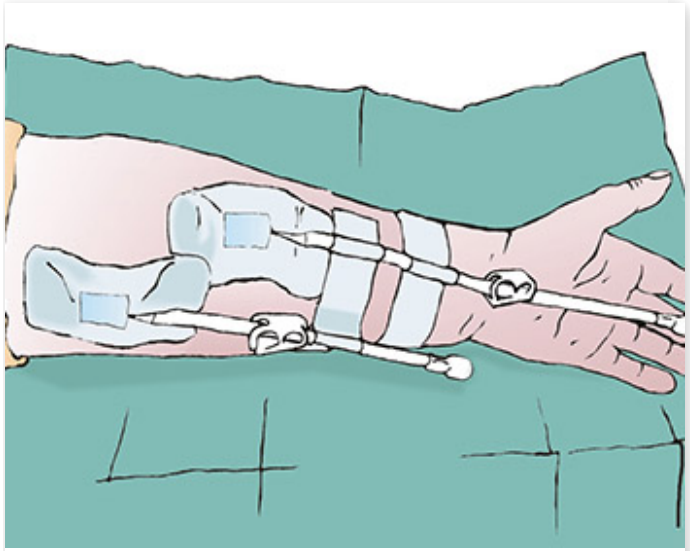

Now the patient is connected to the "artificial kidney".
Dialysis takes 3-4 hours.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026