Matoset dialysis disconnection set – application way

Example of a dialysis disconnection set MA-991-ZEST-488.

1



Before finishing the procedure wash your hands (specialist).

2



Open the dialysis disconnection set.

3



Put one glove on the patient's hand.


©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

4




Then the specialist puts on both gloves.

5

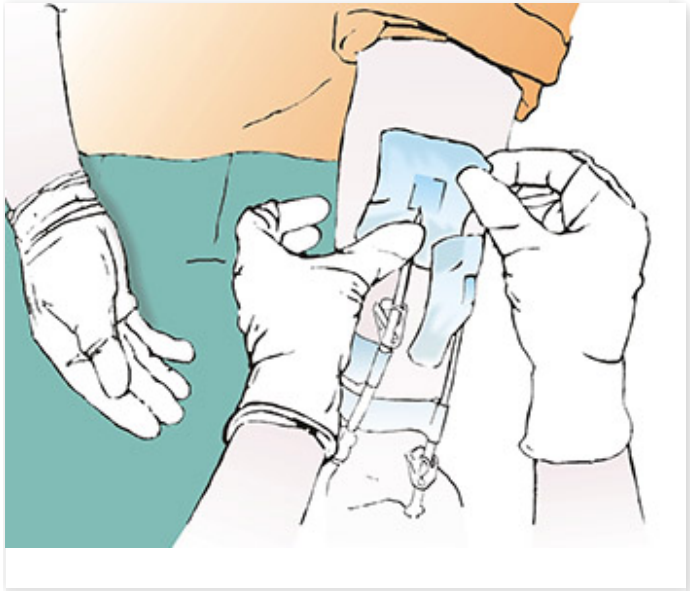

Successively remove the cannulas, using the patient's help to secure the injection sites with the swabs.


6

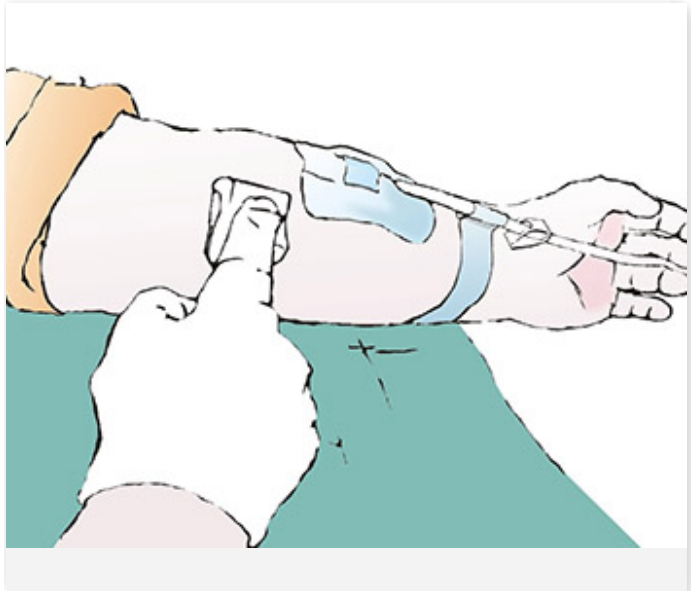

»

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

7

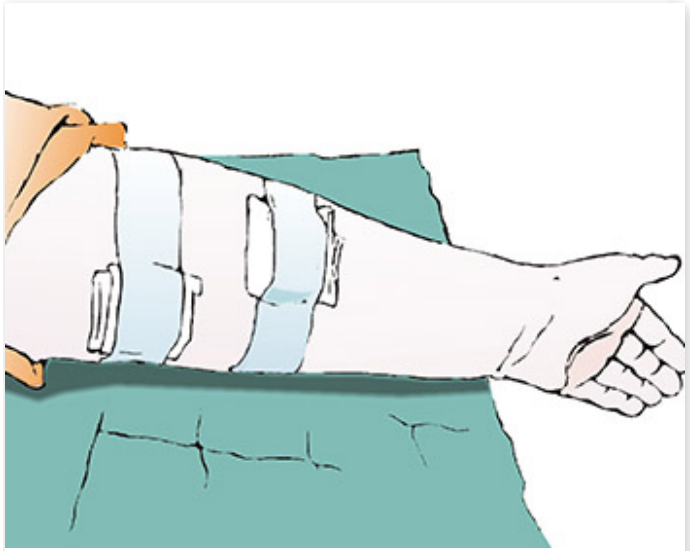

Having removed both cannulas the injection sites must be protected with dressings composed of the nonwoven swabs and adhesive tapes.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026