
Unfolding of Matodrape delivery set

1



Having opened the set put it on the instrumental table. Use the outer drape to secure the table. Put set elements on the table.

2

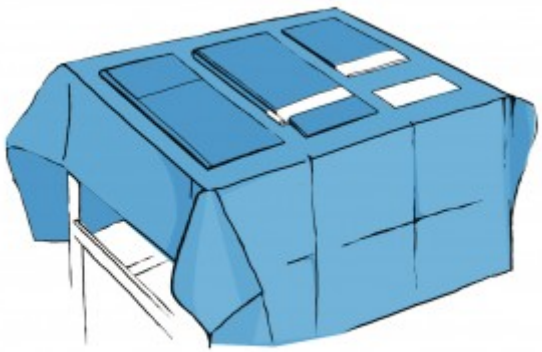

Put the bigger Seni Soft underpad under the patient’s buttocks.


3



Use the drape from the set to cover the patient’s leg to guarantee the sterility of the OP-field.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

4

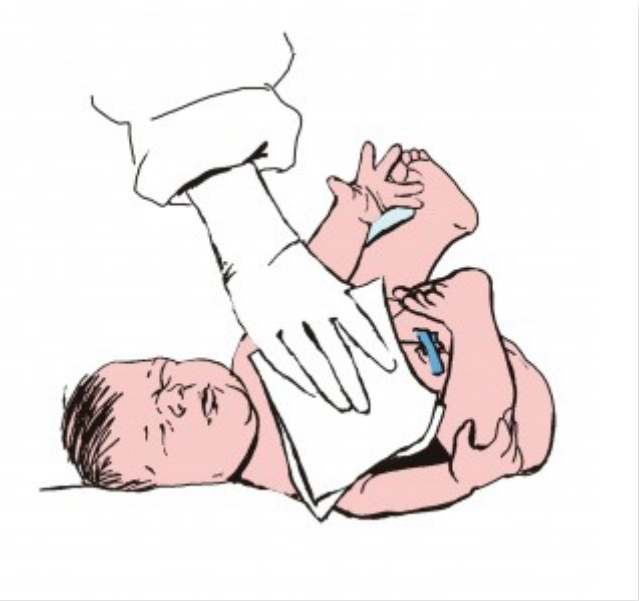

After the delivery dry the baby with white soft nonwoven towels.

5



Cover the baby with smaller Seni Soft underpad and put the ‘baby hat’ on the baby’s head.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026