

# Unfolding of Matodrape laparotomy set

**1**

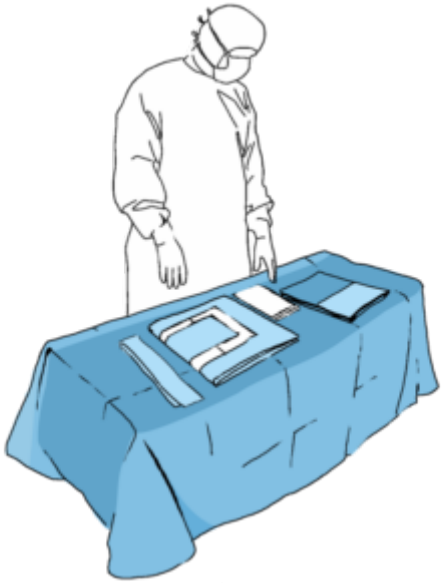

Open the outer packaging in accordance with the principles of asepsis and place the set on the instrument table. Unfold the outer drape on the instrument table and place set elements on it.

**2**

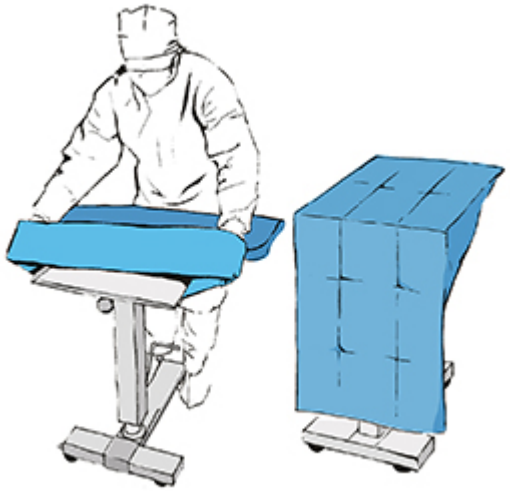

Apply the Mayo table cover.

**3**




Before draping the patient (delimiting the surgical field), pay attention to the direction in which the drape is to be unfolded with the help of the pictograms.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026


# Unfolding of Matodrape laparotomy set

4

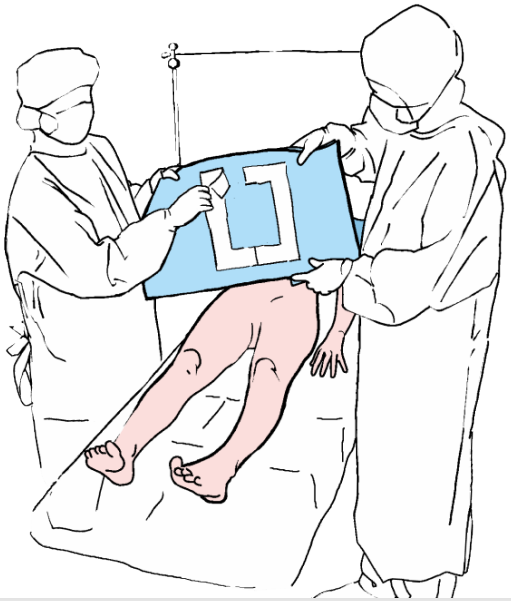


Remove the release paper from the main drape. Adhere the drape to determine the boarders of the surgical field. Then press the adhered area of the drape to increase the adhesion of the adhesive tape to the patient's skin.

5

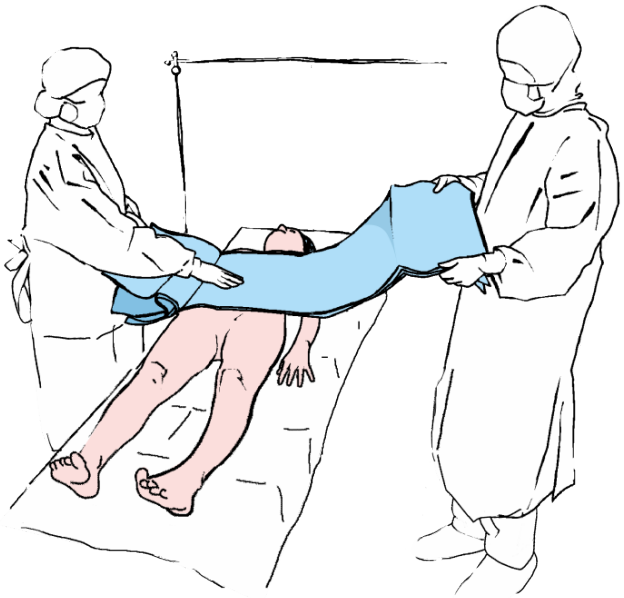

Unfold the sides of the drape so that they hang freely from the operating table.


6

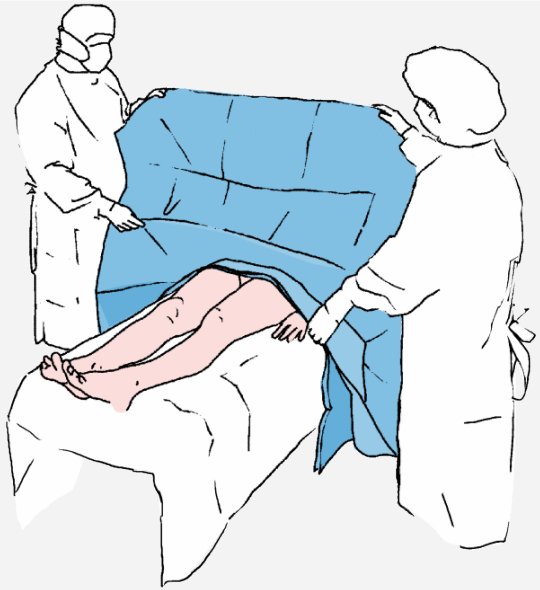

Unfold the top of the drape toward the patient's head (toward the anaesthesia stand). Pass the drape over the frame to create an anaesthesia screen.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026

7

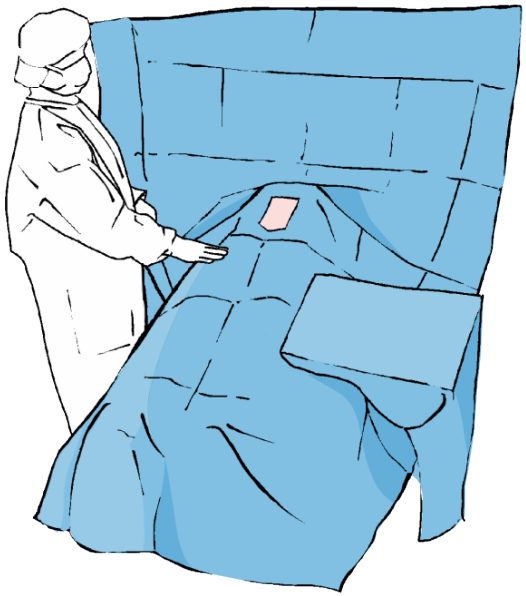

Unfold the lower part of the drape towards the lower end of the operating table to cover the lower limbs of the patient.

©Copyright 2026 en.matopat-global.com | Last updated: 04.11.2026